




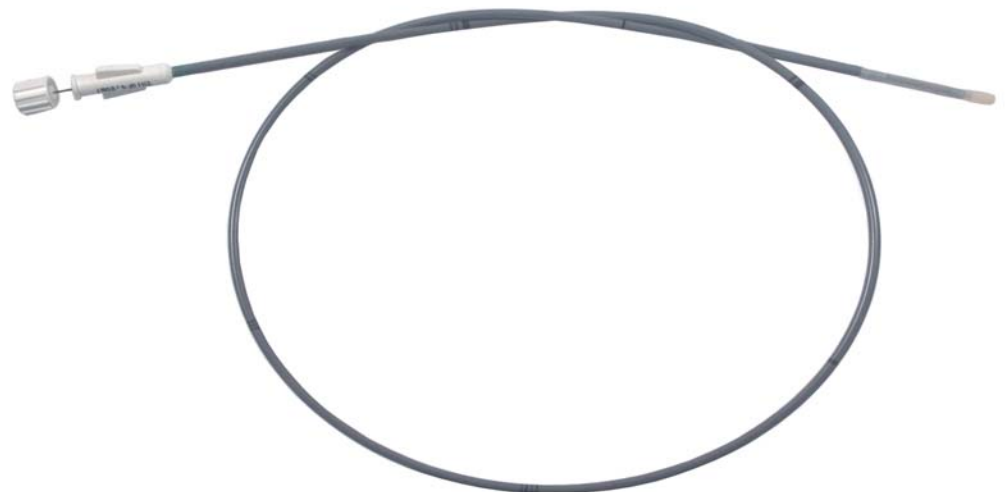


Thrombectomy Catheters - Single Lumen



CODE		COLOUR	EXTERNAL DIAMETER	LENGTH	MAX. INFLATE	BALLOON DIAMETER (MAX.)
310680		Blue	6F (2,00 mm)	80 cm	1,50 ml	12 mm
310880		Brown	8F (2,70 mm)	80 cm	2,25 ml	13 mm
311080		Grey	10F (3,35 mm)	80 cm	4,00 ml	19 mm

Other lengths available upon request



ISOMED Thrombectomy Catheters - single lumen are made with X.R.O. polyurethane. They are colour coded for easy identification. The specially designed latex balloon, fixed with two sutures, is perfectly symmetrical.

Thrombectomy catheters are graduated every 10 cm from the distal end, and are re-inforced with a removable steel wire. The soft and flexible tip of the distal end avoids valve damage.

Double packaged - Sterile (ethylene oxide)