

Jackson-Pratt® Brand Hemaduct® and Hemaduct Gold™ Wound Drains


The Jackson-Pratt® brand Hemaduct® wound drain is uniquely designed with a system of multiple channels and lumens that are interconnected through a series of internal portals. This advanced design offers the following advantages:

- Provides more effective distribution of suction and fluid flow over the entire drain
- Improves collateral drainage by providing alternate pathways around obstructions
- Minimizes aspiration of delicate tissue and clot infiltration with narrow ducts and recessed portals
- Lessens patient discomfort during drain removal by helping prevent tissue ingrowth
- Ensures strength and reliability with unique structural design

Hemaduct drains are made of soft silicone and are available in round and flat versions, with or without a trocar.



Hemaduct round drain

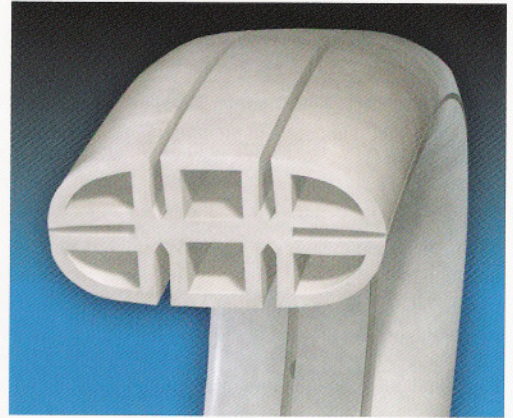

Hemaduct flat drain

Cat. No.	Drain Size*	Ducts	Trocar Size
Hemaduct flat drains without trocar (10/bx)			
JP-HUF073	7mm x 20cm	3/4 (15cm)	
JP-HUF070	7mm x 20cm	Full (20cm)	
JP-HUF103	10mm x 20cm	3/4 (15cm)	
JP-HUF100	10mm x 20cm	Full (20cm)	
Hemaduct flat drains with trocar (10/bx)			
JP-HUF074	7mm x 20cm	3/4 (15cm)	15 Fr x 15cm
JP-HUF071	7mm x 20cm	Full (20cm)	15 Fr x 15cm
JP-HUF104	10mm x 20cm	3/4 (15cm)	15 Fr x 15cm
JP-HUF101	10mm x 20cm	Full (20cm)	15 Fr x 15cm
Hemaduct 3/4-ducted round drains without trocar (10/bx)			
JP-HUR100	10 Fr (1/8") x 40cm	30cm	
JP-HUR150	15 Fr (3/16") x 40cm	30cm	
JP-HUR190	19 Fr (1/4") x 40cm	30cm	
Hemaduct 3/4-ducted round drains with trocar (10/bx)			
JP-HUR101	10 Fr (1/8") x 40cm	30cm	10 Fr x 15cm
JP-HUR151	15 Fr (3/16") x 40cm	30cm	15 Fr x 15cm
JP-HUR195	19 Fr (1/4") x 40cm	30cm	15 Fr x 15cm
Hemaduct fully ducted round drains without trocar (10/bx)			
JP-HUR860	10 Fr (1/8") x 30cm	30cm	
JP-HUR880	15 Fr (3/16") x 30cm	30cm	
JP-HUR900	19 Fr (1/4") x 30cm	30cm	
Hemaduct fully ducted round drains with trocar (10/bx)			
JP-HUR870	10 Fr (1/8") x 30cm	30cm	10 Fr x 15cm
JP-HUR890	15 Fr (3/16") x 30cm	30cm	15 Fr x 15cm
JP-HUR910	19 Fr (1/4") x 30cm	30cm	15 Fr x 15cm

*Length from hub to distal end